

ELECTRIC
MICROMOTOR

DABI ATLANTÉ

ENJOY
THE
STILLNESS

MORE POWER

MORE ACCURACY

LESS NOISE

CONTRA-ANGLE 1:5 LED+

The Contra-Angle operates coupled to the Electric Micromotor, allowing work with maximum rotation of up to 200,000 RPM and torque up to six times higher than a conventional turbine. Via the potentiometer, you are able to select the ideal rotation limit for each type of procedure.

QUADRUPLE SPRAY COOLING

Directed to the drill tip - the double internal spray system, cools the drill cutting area in any working position.

LED+

Optimal illumination through an LED (located on the micromotor) to the contra-angle, via glass fibres. Light power: 25,000 lux

REDUCED HEAD

The Contra-Angle has a small head, with a handpiece angle following the most rigorous international standards - providing greater visibility for the dental professional.



ELECTRIC MICROMOTOR (EMM) CONTRA-ANGLE

DETAILED ACCURACY

INCREASE YOUR VISUAL ACUITY!

A key advantage to the EMM is the LED+ lighting system, diffusing a natural and white light, similar to daylight - superior to traditional LED lighting systems.

Transmission Speed: 1:1
Maximum Speed: 40,000 rpm
Maximum Torque: 3.5 N.cm
Irrigation System: Internal
Lighting system: LED+

Connection Type: ISO 3964
Quick coupling (INTRA)
Weight: 57 grams




ELECTRIC MICROMOTOR BENEFITS

- Work with controlled speeds and high torque.
- Increase your productivity and treatment quality from start to finish - ensuring both you and your patients are satisfied.
- Consistent performance: Work without loss of speed and power, even in high load situations. The Electric Micromotor has constant, reliable operating efficiency.
- Adjust your rotation speed to suit the specific needs of each procedure.
- Preparations with minimal invasiveness and refined details - improving your prosthetic-restorative quality.
- Low noise level: It's electromagnetically induced speed generation, ensures quiet performance.


PORTABLE ELECTRIC MICROMOTOR (TABLE-TOP)

- Operate with total mobility.

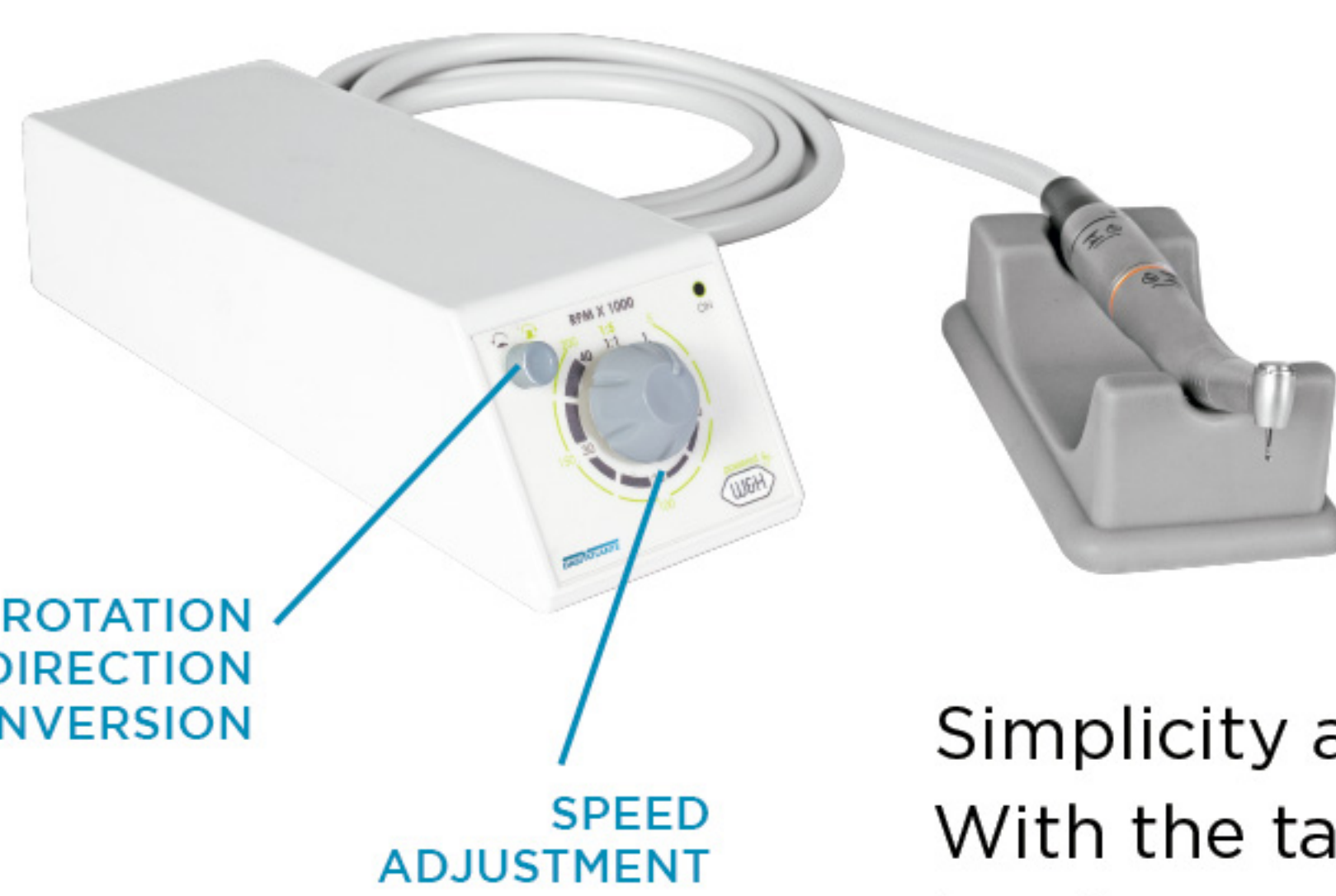

- 1 ON/OFF
- 2 POWER SUPPLY CONNECTION
- 3 BORDEN CONNECTION TERMINAL

Simplicity and practicality:
With the table-top electric micromotor, you can couple the borden terminal of your equipment with ease and use the Electric Micromotor via the equipment's own foot control.

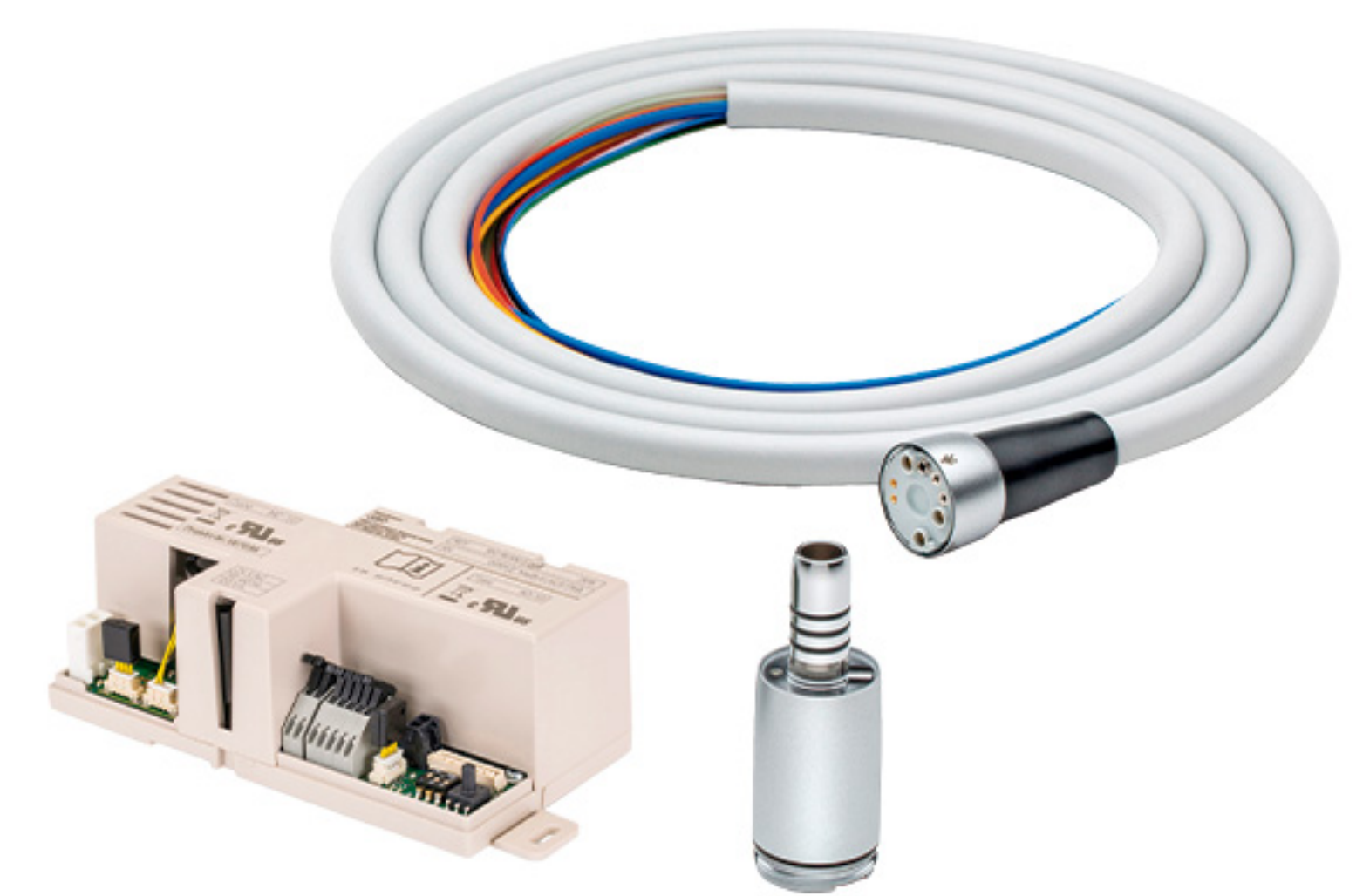

EQUIPMENT KIT ELECTRIC MICROMOTOR (EMM)

The EMM can be purchased as a kit and installed in any dental unit equipment.

THE ELECTRIC MICROMOTOR COMES WITH:

- Electric Micromotor
 - Power Source (table-top EMM)
 - Electronic Module
- *Contra-Angle (Optional)

Illustrative photos. Dabi Atlante reserves the right to make changes to the design and technical specifications of its products, with the authorization of ANVISA, in order to pursue its improvement. Portable Electric Micro Motor: 10101139032. Contra-Angle WH WG-99 LT: 80047309158. The Electric Micromotor Kit is an accessory for the exclusive use of Dental Equipments: 10069210075. 51000000153 - Version May 2022.



dabiatlante



dabiatlante_oficial



dabiatlante

DABI ATLANTE